

The experienced FAMOS and L.R. Instruments partnership brings you a complete range of sterilisation hardware and consumables.



Incorporating both traditional and advanced technologies for the Australian and New Zealand markets.



Build on your experience and take advantage of our complete solution ...

# Consumables for Famos Sterilisation



Total  
Sterilisation  
Confidence  
Full  
Recall  
of Results  
Your Total  
Solution



SASCI Scanning  
SASCI Analysis  
SASCI Storage of  
SASCI Chemical  
SASCI Indicators  
**SASCI**

Daily solution the SASCI will provide you  
and your sterilisation department!



## Helix Bowie Dick System

The efficiency of air removal is a critical parameter for pre-vacuum sterilizers. A steam penetration test should be performed in an empty sterilizer at the beginning of every day to verify and document the efficiency of the vacuum pumps and to indicate possible leaks and/or poor steam quality.

- Get more sensitive reactions with possible air leaks with 6 test pads to indicate pass or failure
- PCD enables life of 1000 cycles (with correct replacement of wear parts)
- Validated according to test pack DIN EN 285

**F15.100.050L** Pack of 250      **F15.100.060L** Pack of 100



## Helix Lot/Batch Control

As the Bowie-Dick Simulation Test is a once-a-day test, performed in an empty autoclave, a lot control test should be used with each batch and documented. The Lot Control Test, by means of a test tube with small lumen, simulates the sterilisation of goods which are difficult to penetrate. Instant results by examining the colour control indicator, means the whole batch can be released immediately after sterilisation.

- Easy to disassemble polypropylene housing which enables quick checking of batch results without waiting to cool

**F15.100.090L** Pack of 100      **F15.100.100L** Pack of 400



## Porous Bowie Dick Test

The benefits of a Traditional Single Use Bowie Dick Test Pack, with the addition of easy and unique Documentation Boxes incorporated on to the sheet and a 5 year shelf life from date of manufacture. The product conforms to EN 867-4 and ISO 11140-4 and is available in packs of 20.

**F15.100.069**



## Class 6 Integrator

Features:

- Conforms to ISO11140-1 Class 6 and EN867 Class D  
(See **ASNZS4187:2003 Part 8.6 Monitoring Sterilizing Cycles**)
- Operates at 121°C & 134°C.
- Incorporates a Clear Reference Image
- Presented in lots of 250 in a practical and user friendly Dispenser System

**F15.100.058**



## Batch/Tracking Sticker

**F15.300.164** Steam  
**F15.300.171** ETO  
**F15.300.205** Hydrogen Peroxide (Plasma)  
**F15.300.172** Formaldehyde  
 Plain Multiple versions available

Double label 37 x 19 mm, 2 times self-adhesive, 750 labels/rolls  
 Class A - EN 867-1



## Batch Labelling Gun for FAMOS Batch Stickers

Very Reliable and incorporate's the following alternatives data

- kind of sterilisation, e.g. steam
- number of sterilizer
- lot number
- personal code for releasing the lot
- date of sterilisation
- expiry date

**F15.300.250**



## Heat Sealers

**F108 SERIES:** An entry level model which is simple, hard wearing, reliable and easy to operate.

**F108, F108T, F108Ti, F108Tx\*** Star Product

**F110 SERIES:** A feature packed medium to high use model suited to hospitals and other large scale health care applications. **F110, F110D, F110DEP, F110DBC**

**F220 SERIES:** The 220 Series heat sealers feature temperature, speed and pressure measuring and monitoring via extra sensors. **F220D, F220DEP, F220DBC**

**Contact Richards Medical for more information on the Famos Heat Sealers.**



## Seal Test

AS/NZS4187 "Informs" the attributes of heat sealing devices. The Seal Test is designed as a daily Evidence Base Quality Assurance check to confirm that the roller pressure and temperature are operating within the correct operation parameters. The result is a quick visual indication of sealer performance and it incorporates simple documentation. Conforms to ISO11607.

**F15.300.038**



## Tape dispensers

**F15.120.026** Metal Housing Single Roll up to 25mm / 50m

**F15.120.027** Metal Housing Two Rolls up to 25mm / 50m

**F15.120.021** Metal Housing Three Rolls up to 25mm / 50m



## Autoclave Tape

**F15.300.006** Standard Tack, 19mm pack of 12

**F15.300.008** Standard Tack, 25mm pack of 9

**F15.300.011** High Tack, 19mm pack of 12

50m/roll with indicator Class A - EN 867, single packed



## Marker Pens

Features:

- Safe to use in the autoclave
- Non-toxic and free from heavy metals
- Instantly touch dry

**CPEN001 BLACK, CPEN002 RED, CPEN003 BLUE, CPEN004 GREEN, CPEN005 ORANGE**



## Steripeel Roll Storage and Cutter

**F50.100.140** TS600 Reel Dispenser & Cutter

**F50.100.409** TS600 Reel Dispenser & Cutter & Storage


**F50.100.080** TS750 Reel Dispenser & Cutter

**F50.100.410** TS750 Reel Dispenser & Cutter & Storage

**F50.100.130** TS900 Reel Dispenser & Cutter

**F50.100.411** TS900 Reel Dispenser & Cutter & Storage

# Typical Sterilisation Process & Quality Assurance Product Options

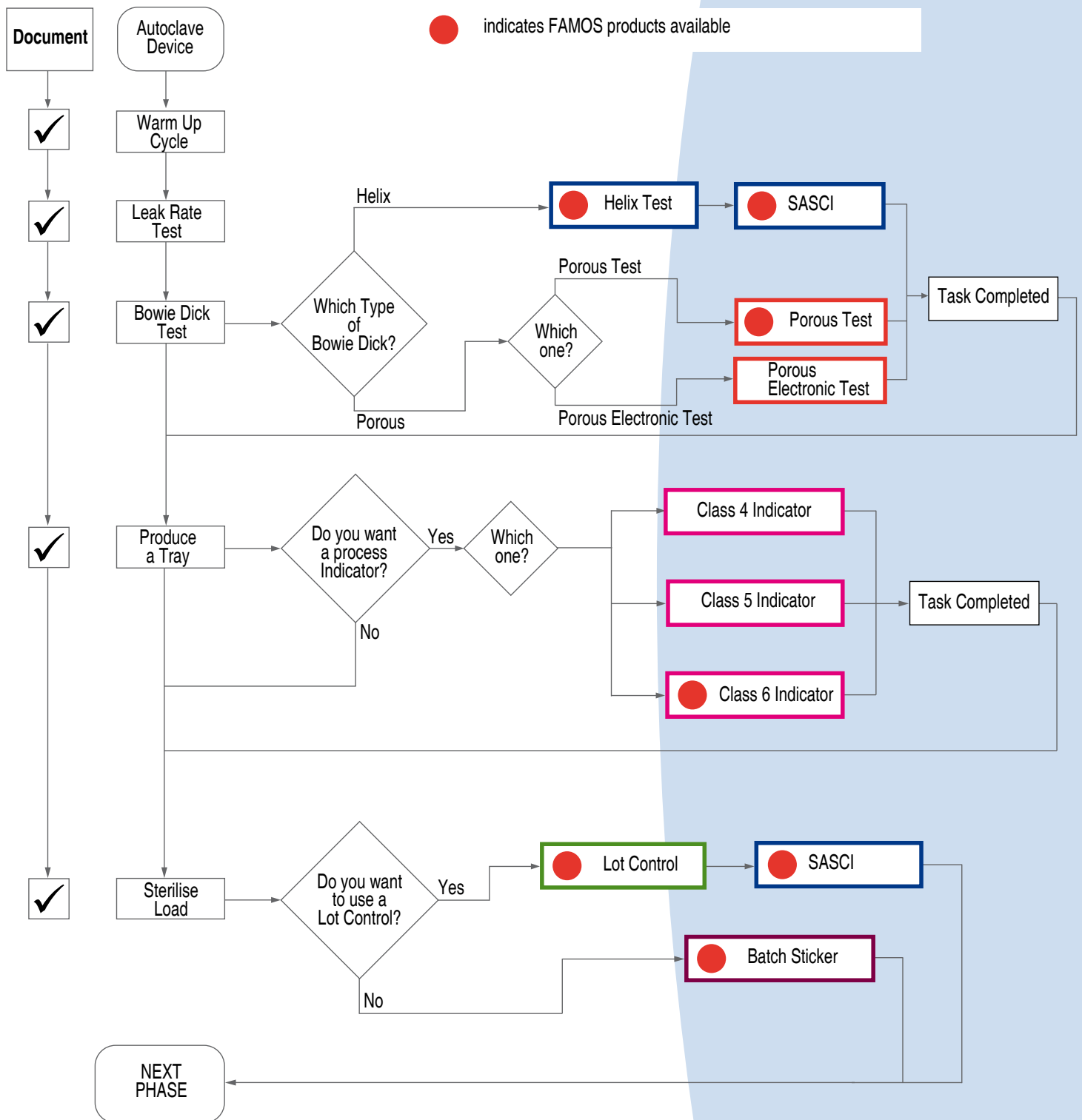

www.famos-med.de



L.R. INSTRUMENTS PTY LTD

UNIT 7, 100 NEW STREET RINGWOOD VICTORIA 3134 AUSTRALIA

Phone: (03) 9870 7066 | Fax: (03) 9870 1638

Email: info@lrinstruments.com.au | Web: www.lrinstruments.com.au

E&OE

First Published: 5 May 2009

Last Modified: 24 November 2009

File Name: brochure\_FAMOS\_CONSUMABLES