

## **Directions for Use**

### **Famos 3 Line Label Gun**

Reference: 15.300.255

#### **Contents**

---


1. General information
2. Setting the print
3. Insert of the label roll
4. Removing the ink roller
5. Inserting a new ink roller

## 1. General information

The Famos 3 Line Label Gun is specially designed to be used with Famos 3 Line Double Labels for documentation and traceability purposes in CSSD.


This Label Gun works very efficiently and it is very easy to use in terms of loading, carrying and printing. Below are the steps also shown in pictures you should follow to make your label gun ready for use.

## 2. Setting the print:

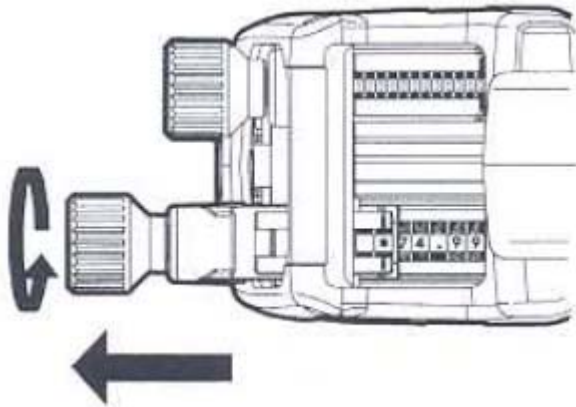


***Pull adjustment knob to the desired position and turn knob to select character to be printed. Each roll has a selection for printing blank.***  
***IMPORTANT: Adjustment knob must be pressed in before operating to avoid damage to the print head.***

## 3. Insert of the label roll:

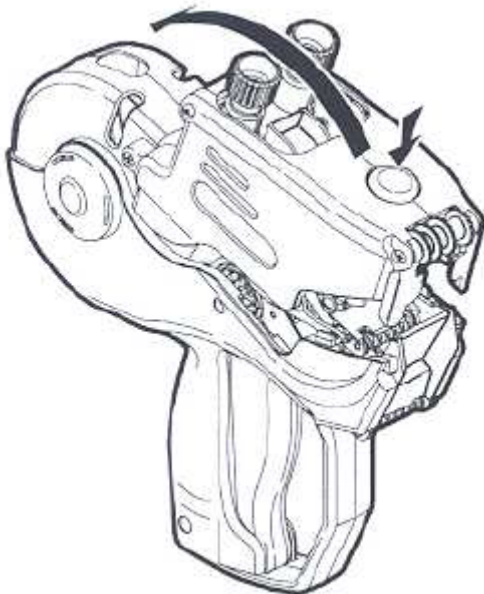


***Push the release button and swing the cover open completely. When you hear a click, don't force the cover anymore.***

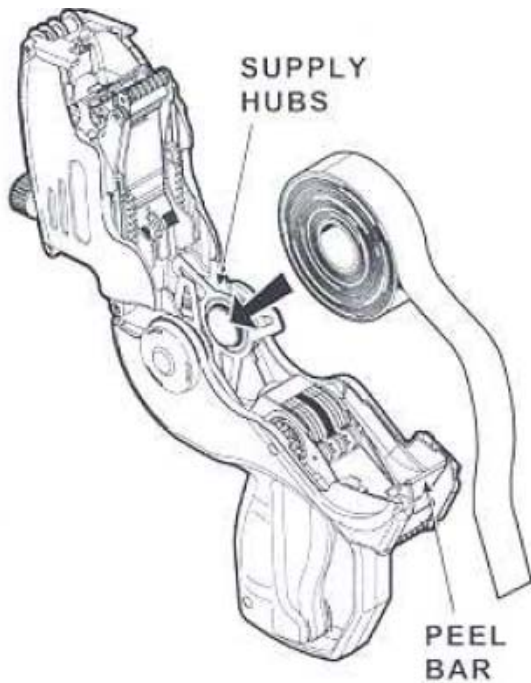

***Please release 12 labels from the roll.***

***Place the label roll between the supply hubs. Drape the supply across the labeller and centre it over the peel bar.***

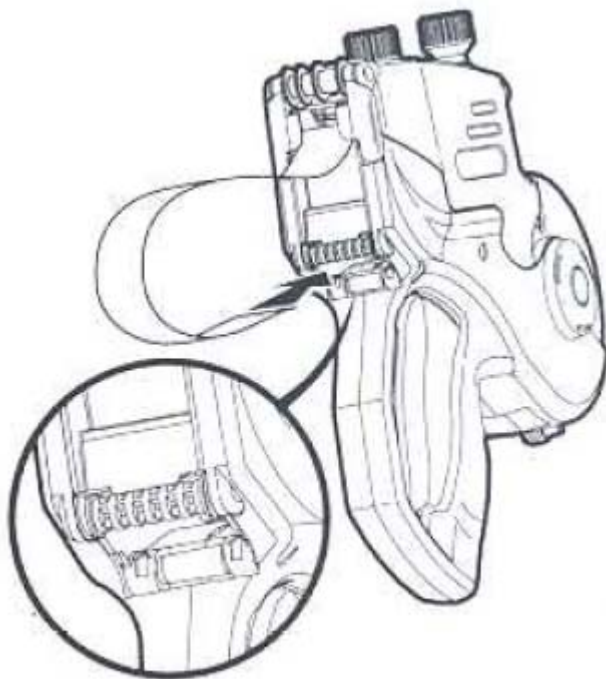


***Close the cover.***

***Feed the label strip under the large roller and into the slot above the arrows.***

***Hold the strip in place and squeeze the trigger until the label strip feed through the rare exit.***

**4. Removing the ink roller:**

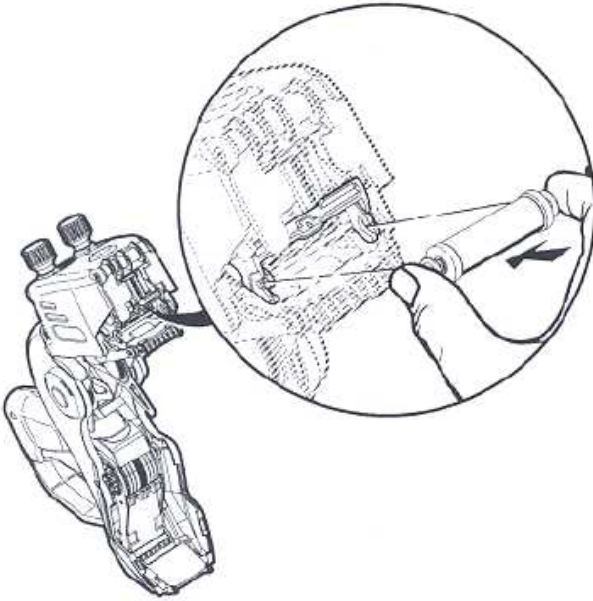


***Open the cover completely.***

***Hold ink roller on each end and push the release bar.***

***Discard the used ink roller.***

**5. Inserting a new ink roller:**

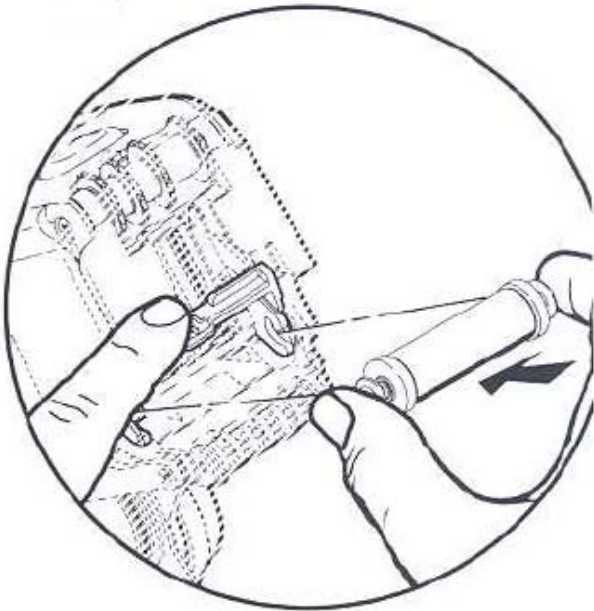

***Open the labeller cover and set aside.***  
***Hold the new ink roller on each end and snap into place.***

***Close the cover.***