


Disposable Manual Resuscitator Adult



Features

- Meets ISO-8382 and ASTM0920 standards.
- Arrow shape texture surface for slip free operation in wet conditions.
- Faster rebound and less memory material applied for ease of operation.
- Ergo designed shape provides ease of handling and less access force.
- Clear resuscitation bag for ease of observation.
- High efficiency silicone duckbill patient valve minimises back leakage.
- Provides low inspiratory and expiratory resistance.
- Swivel 15/22mm patient connector.
- Meets ISO-5356-2 standard.
- 60 cm H₂O pressure relief with lock function for adult model.
- High efficiency reservoir system provides up to 99% oxygen concentration at low flow.
- Swivel elbow oxygen tubing connection for ease of access. Fits-all oxygen connector can be attached to flow meter directly.
- Clear air cushion mask provides better seal, ease of observation during access and also patient comfort.

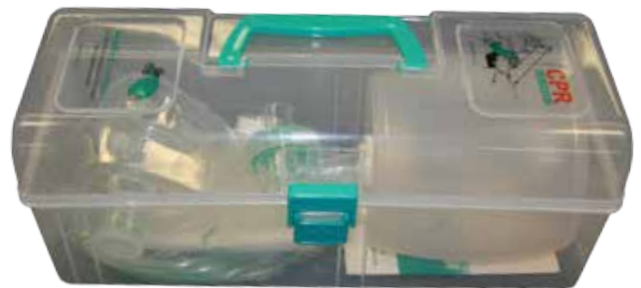

0043



Code	Description	Mask	Reservoir	Pop-off
2268	DISPOSABLE MANUAL RESUSCITATOR ADULT	5305	1000ML	60 CM H ₂ O
0043	SEE-THROUGH CARRY CASE ADULT			



Prepared by: Paragon Medical Pty Ltd trading as LR Instruments ABN 76 133 072 107
 Unit 7/100 New Street (PO Box 458) Ringwood VIC 3134
 Contact Details: P (03) 9870 7066 | F (03) 9870 1638 | E info@lrinstruments.com.au
 File Name: 2268_Disposable-Resuscitator_v3-1

First Prepared: 22/01/15 Eevien Tan
 Release Approval: 22/01/15 Lee Carter
 Version Control: V3.1 Page 1 | 1