



GUEDEL AIRWAYS


Our Guedel Airways are designed for free air flow.

By protecting from airway occlusion from the relaxation of the tongue in an unconscious patient.

Our Guedel Airway devices have a central lumen which allows for various suction catheters to be inserted.

The top section has a colour coded section for rapid identifying of size as well as a thickened area for protection from biting. This bite block is ideal for patients with seizures and unprotected airways. Also it has benefits for the intubated patient disallowing them from crushing the endotracheal tube.

There are various sizes available from neonate to large adult. It is important to choose the right size; the size of the Guedel airway is the length between the corner of the mouth and the ear lobe. It is inserted into the patient's mouth upward while holding the tongue, and rotated when introduced to point downward.

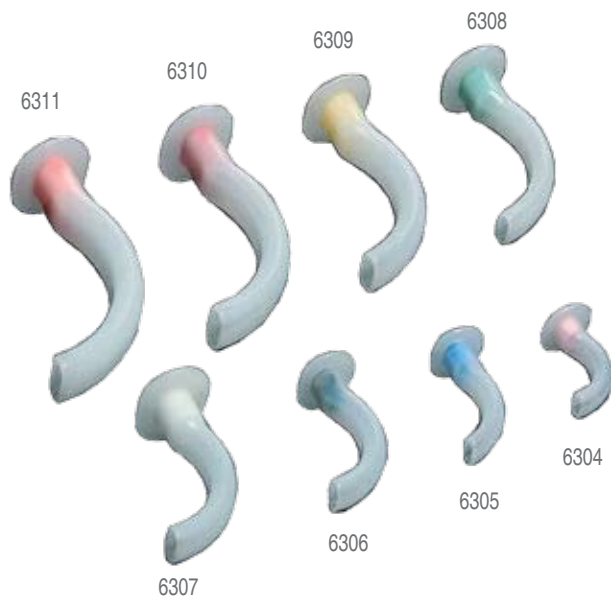

Features

- Features airway and bite block for airway management
- Consists of a large flange, a reinforced bite portion, and a gentle curve that follows the contour of the tongue
- Colour code for size identification conforms to ISO requirements
- Meets ISO-10993:1997 for biocompatibility requirements



CODE	DESCRIPTION
6304	GUEDEL AIRWAY 40 MM COLOUR PINK
6305	GUEDEL AIRWAY 50 MM COLOUR BLUE
6306	GUEDEL AIRWAY 60 MM COLOUR BLACK
6307	GUEDEL AIRWAY 70 MM COLOUR WHITE
6308	GUEDEL AIRWAY 80 MM COLOUR GREEN
6309	GUEDEL AIRWAY 90 MM COLOUR YELLOW
6310	GUEDEL AIRWAY 100 MM COLOUR RED
6311	GUEDEL AIRWAY 110 MM COLOUR ORANGE

E&OE First Published: 20 March 2009
 Last Modified: 24 November 2009
 File Name: info_63XX_GUEDEL-AIRWAYS

